



Intelligent
health systems


CHEMICAL FREE

Mosquito-Click™

If it stings, if it itches...
YOU CLICK!

**DON'T GET CAUGHT OUT
THIS SUMMER!**

A revolutionary remedy
for insect bites!

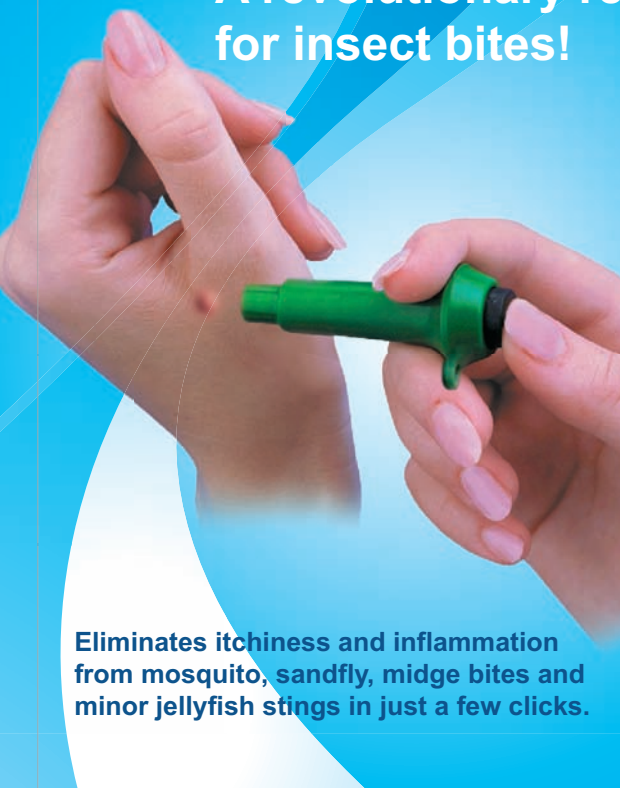

MOSQUITO



SANDFLY
& MIDGE



JELLYFISH



Eliminates itchiness and inflammation
from mosquito, sandfly, midge bites and
minor jellyfish stings in just a few clicks.

UP TO 25,000 CLICKS!

www.mosquitoclick.com



Mosquito-Click™

is a revolutionary and effective method for eliminating irritation and the unsightly effects of insect bites.

Too often we find that our nights are spoilt by mosquitoes. Itchiness resulting from bites persists for a long time and if we give in to temptation to scratch, we end up with skin irritations that can be more serious than the initial bite. To eliminate these effects, a small and very practical device called the Mosquito-Click was invented.

HOW DOES IT WORK?

It functions with piezoelectricity ie. the generation of electric currents by the compression or expansion of quartz and similar crystalline substances. Mosquito-Click therefore does not need a battery to function and its small and practical size makes it convenient to carry along anywhere! Local reactions such as itches, burns, pains and swelling are eliminated in minutes by the reduction or inhibition of histamine. Furthermore, any irritating or poisonous substance introduced under the skin by an insect is altered and therefore inactivated by the application of electric pressure.

INSTRUCTIONS FOR USE

Hold the device like a syringe and position the narrow end on the bite. Press the button a few times. The number of clicks will differ according to the result required.

- 4 to 5 clicks to treat itching
- 10 to 15 clicks to prevent swelling

The number of clicks should be increased for those who are extremely sensitive to mosquito bites. If necessary, the action can be repeated.

The effect of the Mosquito-Click is felt a few minutes after use.

For best results, Mosquito-Click should be used as soon as possible after the bite. It is therefore recommended that you carry it with you, especially during summer. Mosquito-Click is recommended for use on mosquito, sandfly, midge bites and minor jellyfish stings.

WARNING: Do not use together with inflammable substances. Persons wearing cardiac devices or with a heart condition should not use this device. Mosquito-Click can be a choking hazard to children under four, please keep out of reach.

- **CLINICALLY TESTED & CHEMICAL FREE**
- **SUITABLE FOR ADULTS & CHILDREN OVER 12 MONTHS**
- **KEY RING ATTACHMENT**

TECHNICAL SPECIFICATIONS:

- Voltage : 13 kV
- Current : 0.7 mA
- Pulse duration : 0.01 ms
- Diameter : 24mm
- Length : 63 mm
- Weight : 10g
- Power supply : piezoelectricity
- Capacity : 25,000 clicks



Intelligent
health systems

Manufactured by
Tecnimed srl, ITALY for
INTELLIGENT HEALTH SYSTEMS
Factory 2, 8-12 Coora Road,
South Oakleigh Victoria 3166 AUSTRALIA

www.intelligenthealthsystems.com.au

Make an intelligent choice...